


Other materials



Color system

| | |
|-------------------|--------------------------------------|
| Gingiva Composite | Deep Red / Medium Red / Light Red |
| Incisal | Fully Translucent / Semi-Translucent |
| Dentin layer | A1 / A2 / B1 / B2 |
| Opaque | White / Red |

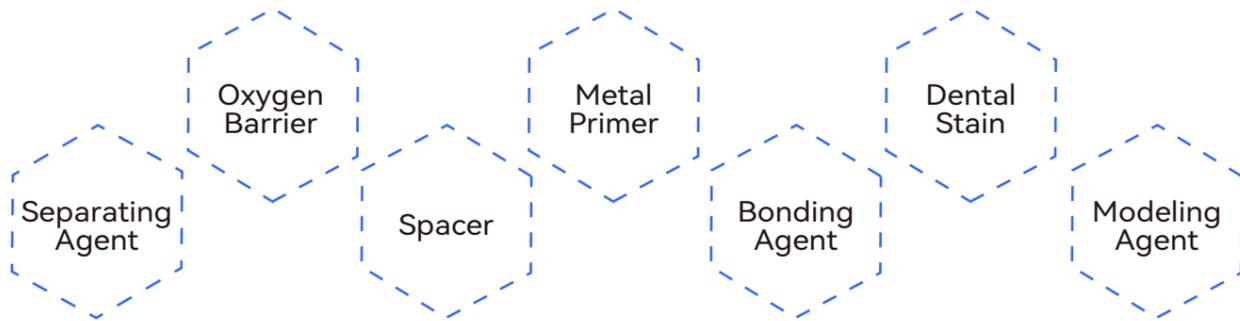

Shenzhen RayForm Technology Co., Ltd.

Nm/Ttl.: Christina(Sales Manager)

Tel.: +86755-27615363 / +8615820462870

E-mail: info@rayformtech.com

Web: www.rayformtech.com

Add.: Building 3 Westside, Guole Science Technology Park, No.1 Lirong Road, Dalang Street, Longhua District, Shenzhen



Ceraform

Dental light-cure composite resin



Why choose us?

Ceraform is high-strength light-curing hybrid ceramic restorative material, combines superior mechanical performance with refined aesthetics. Its optimized resin matrix and finely dispersed ceramic fillers provide excellent strength, wear resistance, and long-term stability under occlusal load. Outstanding polishability, durable gloss retention, natural translucency, and color stability make it suitable for both anterior and posterior restorations, while its excellent handling characteristics ensure precise, efficient, and reliable clinical application.

Indications

Veneer restorations

Repair of ceramic and indirect restorations

Metal- or fiber-reinforced crowns and bridges

Inlays and onlays

Full crowns

Adjustment and repair of conventional acrylic resin teeth

Facing of telescopic crowns and implant-supported restorations



Excellent physical properties

| Physical properties | Ceraform | Ceraform (Flowable) |
|--------------------------------|----------|-----------------------|
| Vicker strength (MPa) | 675 | 385 |
| Compressive strength (MPa) | 362 | 280 |
| Flexural strength (MPa) | 145 | 136 |
| Tooth abrasion (um/200k times) | 0.57 | 0.82 |
| Opaque | 2.2 | 2.5 |

Polymer-ceramic resin series



Dentin layer resin

Flowable capacity: 2g/3g
Paste Capacity: 4g



Gingiva Composite resin

Flowable Capacity: 2g/3g
Paste Capacity: 4g



Translucent resin

Flowable Capacity: 2g/3g
Paste Capacity: 4g

Performance

★ Paste Type ★

- High filler paste formulation enables precise anatomical sculpting and stable build-up.
- Optimized resin system delivers excellent elasticity and toughness, reducing chipping and micro-fractures.
- Superior wear resistance, color stability, and biocompatibility ensure long-term clinical reliability.

★ Flowable Type ★

- Outstanding fluidity allows effortless spreading and uniform thin-layer application.
- Excellent adaptation to complex surfaces supports a wide range of clinical indications.
- High elasticity, wear resistance, color stability, and biocompatibility ensure durable and aesthetic outcomes.



Opaque

Capacity: 2g

Performance

- High masking ability with excellent opacity, effective even in minimal thickness.
- Strong chemical bonding compatibility with light-curing and indirect restorative materials.
- Fast and uniform light-curing performance for stable thin-layer applications.